




Juxta-Fit™

The New Standard
In Comfort and Inelastic Compression


Now, all you have to do is
sit back, relax, and
Juxt-Breathe.

The Juxta-Fit™ Legging System

incorporates state-of-the-art materials with a simple patented design to create one of our most comfortable and effective compression garments ever.

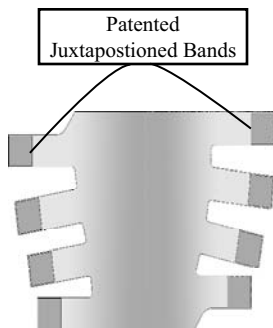


Patented
Juxtapositioned Bands

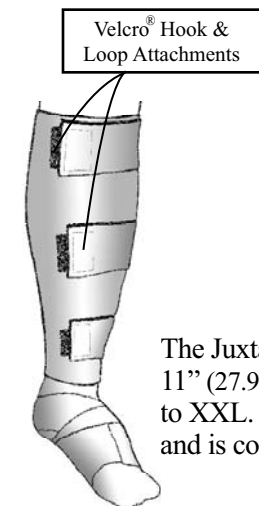

Velcro® Hook &
Loop Attachments

- The Juxta-Lock™ (alternating) Band System allows for quick and easy application and minimal adjustments that help maintain compression levels, a comfortable fit, and therefore increases patient compliance.
- An easy alternative for those who lack the strength or dexterity to use traditional compression stockings.
- The Juxta-Fit™ is made with Breathe-O-Prene™ that allows the garment to breathe keeping the leg cool, dry, and comfortable.
- The material utilizes limited linear stretch technology along the length of the garment that hugs the contours of the leg, reduces slipping, and makes the use of foam inserts virtually unnecessary.
- Adjusting the bands of the garment to a firm and comfortable compression achieves therapeutic levels of gradient compression and can be worn 24 hours* a day with the bands being slightly loosened for night time wear.
- The Juxta-Fit™ contains SilverTec™ in the Breathe-O-Prene™ liner, which is an anti-microbial, anti-static, thermodynamic, anti-odor agent that promotes healthier skin and faster healing.
- The Juxta-Fit™ is specially designed to use varying bandwidths to control and assure gradient compression throughout the entire length of the garment regardless of limb shape.

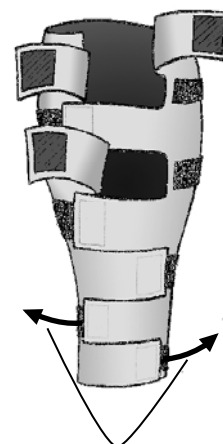

Pulling Bands with
Equal Tension Assures
Gradient Compression

The Standard

The Juxta-Fit™ Lower Legging is available in two standard lengths, 11" (27.94 cm) and 14" (35.56 cm), and in sizes ranging from Small to XXL. It can also be ordered custom to fit almost any leg shape and is covered by a one year limited warranty.

The Custom

The Juxta-Fit™ Whole Legging, Upper Legging, and Knee Piece are custom made only, and can be combined to create six different garment options. Whether you need a garment that covers from just above the ankle to above-the-knee, or to cover from below the knee up to the groin, the Juxta-Fit™ can do it. They are also all covered by a limited one year warranty.

Each custom made garment can be designed to fit the shape and contours of almost any leg and is built to accommodate a 20% circumference reduction in the limb and **still fit comfortably.**

*The Juxta-Fit™ Whole Legging, Knee Piece, Lower Legging w/Knee, Upper Legging, and Upper Legging w/Knee are recommended for night-time use only.

A Physician's prescription is required.

Go to our website:
www.circaid.com
for full color pictures
Phone: (800) CIRCAID (247-2243)

